



## SeroELISA™ Chlamydia IgG

Enzyme-Linked Immunosorbent Assay (ELISA) for the determination of specific IgG antibodies to *Chlamydia* in human serum

### Instruction Manual

Test kit for 96 determinations  
(Catalog No. 111-01)

For In Vitro Diagnostic Use  
For professional use only  
Store at 2-8°C. **Do Not Freeze**

### Savyon® Diagnostics Ltd.

3 Habosem St. Ashdod 7761003  
ISRAEL

Tel: +972.8.8562920

Fax: +972.8.8523176

E-mail: [support@savyondiagnosics.com](mailto:support@savyondiagnosics.com)

### Intended Use

The SeroELISA™ Chlamydia IgG kit is intended for the determination of specific IgG antibody to Chlamydia in a single human serum sample or to evaluate paired sera, by an Enzyme-Linked Immunosorbent Assay (ELISA).

### For In Vitro Diagnostic Use.

### Introduction

Chlamydia is a gram-negative obligate intracellular bacteria that causes acute and chronic diseases in mammalian and avian species. The genus Chlamydia is comprised of four species: *C.trachomatis*, *C.pneumoniae*, *C.psittaci* and *C. pecorum*.

*C.trachomatis* is divided into 15 serovars (1). Serovars A, B, Ba and C are agents of trachoma (2), the leading cause of preventable blindness endemic in third world countries. Serovars L<sub>1</sub>-L<sub>3</sub> are the agents of lymphogranuloma venereum. Serovars D-K are the common cause of sexually transmitted genital infection worldwide: cervicitis, endometritis/salpingitis (3) in females and urethritis (4) in both males and females. Endometritis/salpingitis can lead to tubal occlusion with a higher risk of extra uterine pregnancy and infertility. Genital infection may cause an acute and persistent infection occasionally without any clinical symptoms. Generally, these infections are treatable once they are diagnosed. However without any treatment the infection may progress to a severe chronic inflammation leading to infertility, ectopic pregnancy, induced abortion or premature delivery. Moreover, infants to infected

mothers may be infected during birth, leading to conjunctivitis or pneumonia (5).

*C.pneumoniae* is an important respiratory pathogen in humans and causes up to 10% of community-acquired pneumonia cases. It has been associated with acute respiratory diseases, pneumonia, asthma, bronchitis, pharyngitis, acute chest syndrome of sickle cell disease, coronary heart disease, and Guillain-Barre syndrome (6).

*C.psittaci* infects a diverse range of host species from molluscs to birds to mammals and also causes severe pneumonia.

Serodiagnostic tests, which rely on specific immunologic markers, serve as a non-invasive diagnostic tool in identification of both distal and deep infections (7). It has been found that elevated titers of anti-Chlamydia IgG antibody are indicative of active chlamydial infection (3,8).

High levels of anti-Chlamydia IgG antibody are of diagnostic value in chronic or systemic infections such as salpingitis, mechanical infertility, perihepatitis, epididymitis, Reiter's syndrome and pneumonitis (3, 9, 10).

Savyon Diagnostics has developed an enzyme-linked immunosorbent assay for the detection of specific IgG antibody to Chlamydia.

SeroELISA™ Chlamydia test employs the L<sub>2</sub> serovar broadly reacting antigen of *C. trachomatis*. It will detect *C. trachomatis*, *C. psittaci* and *C. pneumoniae* (TWAR) antibodies.

### Principle of Test

- Human serum to be tested is brought into contact with the antigenic material coating the microtiter wells. Specific antibody, if present in the patient serum will bind to the attached antigen, a complex is formed and all the serum components are washed away in the wash phase.
- Horseradish peroxidase (HRP) conjugated anti-human IgG (γ chain specific) is added to the wells. If an antigen-antibody complex was formed in the previous step, the peroxidase conjugated antibody will bind to the antibody moiety of the complex. If no antigen-antibody complex was formed in the previous step, the conjugate is washed away in the wash phase.
- TMB-Substrate is added. A positive reaction is indicated by a blue to deep blue color which develops in the test wells following enzymatic reaction of the peroxidase moiety with peroxide and the chromogen reactant. After the enzymatic reaction is stopped by an acidic solution. The absorbance of the test wells is determined at 450nm by a spectrophotometer.
- The absorbance at 450nm is indicative of IgG anti-Chlamydia titer in patient serum specimens.

### Assay Procedure

1. Chlamydia Antigen Attached to Solid Phase (Ag) + Serum positive for IgG Anti-Chlamydia (Ab<sub>1</sub>)  
↓  
AgAb<sub>1</sub> Complex
2. AgAb<sub>1</sub> Complex + HRP Conjugated Anti-Human IgG (Ab<sub>2</sub>)  
↓  
AgAb<sub>1</sub>Ab<sub>2</sub> Complex
3. AgAb<sub>1</sub>Ab<sub>2</sub> Complex + TMB-Substrate  
↓  
Blue Solution  
↓ ← Chromogen Stop Solution  
Yellow Solution  
(Absorbance Determination at 450nm)

### Warning and Precautions

- **Warning:** THE CHLAMYDIAL ANTIGENIC MATERIAL HAS BEEN INACTIVATED AND CONTAINS NO DETECTABLE LIVE ORGANISMS. HOWEVER, THE STRIPS SHOULD BE HANDLED AND DISPOSED OF AS WOULD ANY POTENTIALLY BIOHAZARDOUS LABORATORY MATERIAL.  
**Precautions:** This kit contains human sera which have been tested by FDA approved techniques, and found to be negative for HBV antigen and for HCV and HIV 1 and HIV 2 antibodies. Since no known method can offer complete assurance that products derived from human blood do not transmit infection, all human blood components supplied in this kit must be handled as potentially infectious serum or blood according to the recommendations published in the CDC/NIH manual "Biosafety in Micro Biological and Biomedical Laboratories, 1988".
- Substrate/Chromogen Solution is an irritant material to skin and mucous membranes. Avoid direct contact.
- **For *In-vitro* Diagnostic Use.**

### Kit Contents

1. Precoated microtiter plate (96 wells per frame). Each sachet contains one microtiter plate comprising 12 removable strips in a plastic frame. Each strip is coated with Chlamydia antigens.  
**1 Unit**
2. Positive Control (human serum positive for anti-Chlamydia IgG antibody). Ready-to-Use.  
**1 Vial, 2.0ml**
3. Low Positive Control (human serum low positive for anti-Chlamydia IgG antibody). Ready-to-Use.  
**1 Vial, 2.0ml**

4. Negative Control (human serum negative for anti-Chlamydia IgG antibody). Ready-to-Use.  
**1 Vial, 2.0ml**
5. Ready To Use HRP Conjugated Anti-Human IgG ( $\gamma$ -chain specific).  
**1 Vial, 10ml**
6. Serum Diluent, Ready-to-Use.  
**2 Bottle, 60ml**
7. Concentrated Wash Buffer (x20).  
**1 Bottle, 100ml**
8. TMB Substrate, Ready-to-Use.  
**1 Vial, 16ml**
9. Stop Solution (1M H<sub>2</sub>SO<sub>4</sub>), Ready-to-Use.  
**1 Bottle, 16ml**
10. Plate cover. **1**
11. Instruction Manual. **1**

### Materials Required But Not Supplied

1. Clean test tubes for dilution of patients sera.
2. Adjustable micropipettes, or multichannel pipettes (5-50, 50-200 and 200-1000 $\mu$ l ranges) and disposable tips.
3. Disposable plastic pipettes (assorted sizes) and safety pipetting devices.
4. One liter volumetric flask.
5. One 50ml volumetric cylinder.
6. ELISA plate washer or wash bottle.
7. Paper towels or absorbent paper.
8. Vortex mixer.
9. A 37°C water bath with a lid, or a moisture chamber placed in a 37° ± 1°C incubator.
10. ELISA-reader with 450nm filter.
11. Distilled or double deionized water for the dilution of Concentrated Wash Buffer.

### Storage and Shelf Life of Reagents

All the materials supplied should be stored at 2° to 8°C. If kept at 2° to 8°C the test reagents are stable until the expiry date indicated on the kit pack. Exposure of originally stoppered or sealed components to ambient temperature for a few hours will not cause damage to the reagents. **DO NOT FREEZE!**

When a kit is in use, the shelf life of the original material is 60 days from the day first opened. Once opened, the aluminum foil sachet containing the strips should be sealed with a tape. The dehydrating packet should not be removed.

### Specimen Collection

Serum specimens should be collected aseptically and stored at 2° to 8°C with 0.05% sodium azide (NaN<sub>3</sub>) as a preservative if they are to be tested within a few days. For longer periods, aliquots of serum specimens should be stored at -20°C. Since turbid or hemolytic serum samples may give less reproducible results, it is highly recommended that serum samples tested be clear and non hemolytic.

## Assay Procedure

### Notes:

- a) The components of this kit have been tested as one unit. Do not mix components from different kit lots or other manufacturer's kits.
- b) All reagents should reach room temperature before use. Serum Diluent and Conjugate Diluent gelatinize when refrigerated. If needed, accelerate liquefaction by warming these components at 37°C for several minutes. Salt crystals may form in the Concentrated Wash Buffer when stored at 2° to 8°C. These crystals should be completely redissolved by warming at 37°C before dilution.
- c) Do not perform the test in the presence of reactive vapors (e.g. from acid, alkaline or aldehyde substances) or dust, since the enzymatic activity of the HRP conjugated Anti-Human IgG may be affected.
- d) Do not touch the top of the strips. Do not touch the edges of the wells with the pipette tips when dispensing reagents.
- e) Use disposable pipette tips. Avoid cross contamination between reagents.
- f) Tap vial lightly on hard surface to free liquid that might be entrapped in the cap.
- g) Avoid entrapping air bubbles in the wells.
- h) Dispense liquids slowly to avoid spraying.
- i) Positive Control, Low Positive Control and Negative Control sera should be run together with serum specimens each time the test is performed.
- j) One well per test should be used for a blank value each time the test is performed.
- k) All the procedure steps should be performed sequentially without interruption.

### Test Procedure -Manual

#### Automation protocol available upon request

#### A) Washing of Strips

*Prewashing of strips is advisable but not mandatory. In any event, strips should be wet by Wash Buffer solution before application of test samples and dried by tapping them over clean absorbent paper prior to testing.*

*If automatic ELISA plate washer is unavailable procedure should be performed as follows:*

1. Remove required number of strips from their aluminum foil sachets and insert them in the plate frame.
2. Dilute the Concentrated Wash buffer 1:20 with distilled water.  
*For Example:* For one strip prepare 100ml of Wash Buffer (5ml Concentrated Wash Buffer with 95ml distilled water). Swirl gently for 20 minutes. Wash Buffer working solution should be prepared before use and excess discarded.
3. After incubation fill each well with Wash Buffer solution up to the middle of the well.
4. Allow to soak for 2 minutes, then discard the strip contents. Repeat these steps **five** times.
5. Dry the top of the strips and frame by gently tapping them over clean absorbent paper.

*Complete washing of wells after incubation is essential for best results. Traces of Wash Buffer should not be left in the wells.*

#### B) Incubation of Serum Samples and Controls

6. Dilute each patient serum 1:256 with Serum Diluent as follows:

1:16        -Add 10µl of patient serum to 150µl of Serum Diluent

1:256       -Add 10µl of 1:16 dilution to 150µl of Serum Diluent

**Controls** are provided in a ready-to-use form and should not be diluted.

7. Pipette 50µl of ready-to-use Positive control, Low Positive Control, Negative Control and 1:256 patient serum dilution into corresponding separate wells of the test strips. Pipette 50µl of Serum Diluent into one well for blank value.

*Pipetting of controls and serum specimens into the wells should not exceed 10 minutes.*

8. Cover the strips with a plate cover and incubate for 30 minutes at 37°C in a moisture chamber.
9. Discard the liquid content of the wells. Wash the wells **five** times and dry as in steps A) 3-5.

#### C) Incubation with Conjugate

10. Pipette 50µl of ready to use HRP conjugated Anti-Human IgG into each well.
11. Cover the strips with a plate cover and incubate for 30 minutes at 37°C in a moisture chamber.
12. Discard the liquid content of the wells, wash them **five** times and dry as in steps A) 3-5.

#### D) Incubation with TMB-Substrate

13. Pipette 100µl of TMB-Substrate into each well.
14. Cover the strips with a plate cover and incubate at room temperature for 15 minutes.
15. Stop the reaction by adding 100µl chromogen stop solution to each well. Pipette the chromogen stop solution in the same sequence and at the same time intervals as the TMB Substrate in step D14.
16. Calibrate the spectrophotometer on the blank well. Determine the absorbance at 450nm and record the results.

*Immediate determination of the absorbance is advisable but not mandatory. Absorbance determination should not exceed 30 minutes, following stopping of chromogenic reaction.*

## Acceptability of Criteria of the Test

A test-run is valid if:

- Positive Control absorbance is  $\geq 0.8$  at 450nm
- Low Positive Control absorbance is  $\geq 0.4$  at 450nm
- Negative Control absorbance is  $\leq 0.15$  at 450nm

If these conditions are not fulfilled, the test run is invalid and should be repeated.

### Calculation of Cut-Off Value (COV)

The cut-off value is calculated according to the following formula:

$$\text{COV} = 0.154 \times (\text{Pc} - \text{Nc}) + \text{Nc}$$

Pc = Absorbance of Positive Control at 450nm

Nc = Absorbance of Negative Control at 450nm

### Interpretation of Test Results

Absorbance at 450nm	Interpretation of Results	Estimated Chlamydia IgG Titer
Below COV – 0.03	Negative	<256
COV ± 0.03	Equivocal*	256
Above COV +0.03 up to Low Positive Control	Low Positive	256-512
Above Low Positive control up to Positive Control	Positive	1024-2048
Above Positive Control	High Positive	≥2048

\* Retest serum samples classified equivocal. If equivocal result is repeated testing of subsequent serum sample is recommended.

### Significance of Results

Clinical experience and recent studies indicate that high levels of IgG antibody to Chlamydia may serve as a marker for active or chronic Chlamydia infection.

### Endpoint Titer Determination

For monitoring of the immune profile of a patient, acute stage and convalescent serum samples should be compared. Endpoint titration of paired sera should be performed with both serum samples tested in the same run. The highest serum dilution with absorbance above the test cut-off value is the **end titer**.

### Limitations of Assay

- No single serological test should be used as the only criterion for diagnosis. All clinical and laboratory data should be taken into account.
- The test is a single serovar (L<sub>2</sub>) ELISA. L<sub>2</sub> contains antigenic determinants existing in serovars of *Chlamydia trachomatis* as well as the group antigen. Antibodies against *Chlamydia psittaci*, *Chlamydia pneumoniae* (TWAR) and *Acinetobacter calcoaceticus* may be detected by this ELISA.

- The significance of specimen titers must be considered in relation to the characteristics of the population being tested. These characteristics would include age, geographical location and sexual behavior, among other factors.
- This test will not indicate the site of chlamydia infection(s). It is not intended to replace cell culture isolation, if available.
- The absence of measurable serum antibodies does not exclude the possibility of a chlamydial infection.
- Bacterially contaminated or hyperlipaemic serum may cause erroneous results.

### Performance Characteristics

SeroELISA™ Chlamydia test was compared to IPAzyme™ Chlamydia IgG/IgA test (Savyon Diagnostics product, Cat. No. 011-01) which is an accepted serological test for detection of IgG antibody to Chlamydia.

The study population included patients with suspected Chlamydia infections as well as healthy individuals (n=328). The results are summarized as follows:

#### Comparison of SeroELISA™ with IPAzyme™

SeroELISA™ \ IPAzyme™	Positive	Negative	Total
Positive	116	8	124
Negative	10	194	204
Total	126	202	328


Overall agreement: (310/328) x 100 = 94.5%

### Cross Reaction

Hospitalized patients, infected with *Neisseria gonorrhoea*, *Staphylococcus aureus* and *Peptostreptococcus anaerobius*, who were diagnosed by commercial serology kits, were also tested with the SeroELISA Chlamydia kit. There was no significant cross-reaction detected.

## Bibliography

1. **Yuan, Y., Zhang, Y. X., Watkins, N. G. and Caldwell, H.D.** (1989). Nucleotide and Deduced Amino Acid Sequences for the Four Variable Domains of the Major Outer Membrane Proteins of the 15 *Chlamydia trachomatis* Serovars. *Infection and Immunity*. **57**:1040-1049. Copyright 1989, American Society for Microbiology.
2. **Trehanne J. D.** (1985). The community epidemiology of trachoma. *Rev Infect Dis*. **7**:760-763.
3. **Piura, B., Sarov, I., Sarov, B., Kleinman, D., Chaim, W. and Insler, V.** (1985). Serum IgG and IgA antibodies specific for *Chlamydia trachomatis* in salpingitis patients as determined by the immunoperoxidase assay. *Eur. J. Epidemiol* **1**: 110-116.
4. **Wang, S.P., Grayston, J.T., Kuo, C.C., Alexander, E.R., and Holmes, K.K.** (1977). SeroDiagnosis of *Chlamydia trachomatis* infection with the microimmunofluorescence test. In: *Nongonococcal urethritis and related infection*, D. Hobson and K.K. Holmes (Eds), P. American Society for Microbiology, Washington DC. p. 237-248.
5. **Thompson III S. E., and Dretler R. H.** (1982). Epidemiology and Treatment of Chlamydial Infections in Pregnant Women and Infants. *Review of Infectious Diseases* **4**:S747
6. **Saikku P., Mattila, K., Nieminen, M.S., Huttunen, J.K., Leinonen, M., Ekman, M.R., Makela, P.H., and Valtonen, V.** (1988). Serological evidence of an association of a novel Chlamydia, TWAR, with chronic coronary heart disease and acute myocardial infection. *Lancet* **II**: 983-986.
7. **Sarov, I., Kleinman, D., Cevenini R., Holcher, G., Potashnik, G. Sarov, B. and Insler, V.** (1986). Specific IgG and IgA Antibodies to *Chlamydia trachomatis* in Infertile Women. *Int. J. Fertil* **31**:193-197.
8. **Csango, P.A., Sarov, B., Schiotz, H., Sarov, I.** (1988) Comparison Between Cell Culture and Serology for Detecting *Chlamydia trachomatis* in women seeking abortion. *J. Clin. Pathol.* **41**:89-92.
9. **Tsunekawa, t. and Kumamoto, Y.** (1989). A Study of IgA and IgG Titers of *C.trachomatis* in Serum and Prostatic Secretion in Chronic Prostatitis. *J. J. A. Inf. Dis.* **63**:130-137.
10. **Kletter, Y., Caspi. D., Yarom, M., Sarov, B., Sarov, I. And Tanay, A.** (1988) Serum IgA and IgG Antibodies Specific to *Chlamydia* in Patient's with Reiter's Syndrome (RS). In: *Proceedings of the European Society for Chlamydia Research, Societa Editrice Esculapio, Bologna.* p. 170.

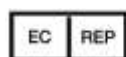

**Obelis** s.a. (European Authorized Representative Center)  
Boulevard Général Wahis 53, 1030 Brussels, Belgium  
Tel.: +32.2.732.59.54 Fax.: +32.2.732.60.03  
e-mail: mail@obelis.net